

# KaVo ProXam iCam

New intraoral camera

# KaVo ProXam iCam

- Intraoral camera with high image quality
- Plug & Play (USB camera)
- Lightweight and convenient size
- Treatment unit integration (1,3m cable) or stand-alone use (2,5m cable)
- Easy operation via camera buttons or KaVo foot control
- Hygienic sheaths to protect camera from contamination



# KaVo ProXam iCam

## Intraoral Mode


- High image quality and depth of focus
- Macro mode for close-up images



# KaVo ProXam iCam

## Operation

- Two camera buttons for intuitive use
- Status LED
- KaVo foot control integration:
  - KaVo treatment unit foot control
  - Stand-alone foot control

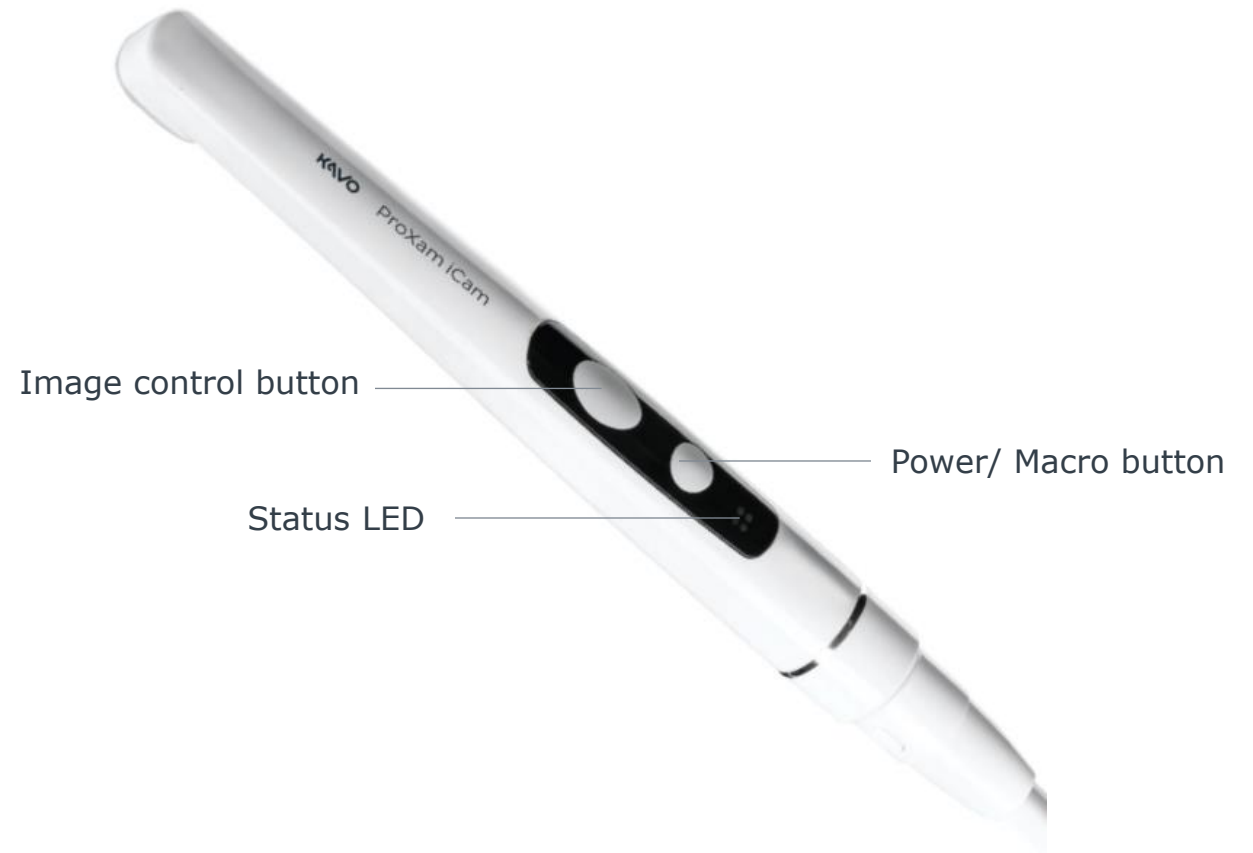

# KaVo ProXam iCam

## Configuration options and accessories

- Two versions of the camera available:
  - KaVo ProXam iCam 1,3m for treatment unit integration (Mat.no. 2.023.8000)
  - KaVo ProXam iCam 2,5m for stand-alone use (Mat.no. 2.023.8100)
- Wall mount included in standard scope of delivery

### **Accessories:**

- Camera Sheaths, 500 pcs. to cover the camera handpiece (Mat.no. 2.023.8002)
- Foot Control Vision Full HD (Mat.no. 1.013.4952)
- Adapter handpiece standing for better placement in KaVo treatment units (Mat.no. 1.011.0275)
- Cradle Vision Full HD for better placement e.g., on a table (Mat.no. 1.013.5500)

# KaVo ProXam iCam

## Software integration

Seamless integration in KaVo patient communication systems:



| Treatment unit integration                                                                       |                                            | Stand-alone use                                                          |                                                     |
|--------------------------------------------------------------------------------------------------|--------------------------------------------|--------------------------------------------------------------------------|-----------------------------------------------------|
| KaVo CONNECTbase                                                                                 | KaVo CONEXIOcom                            | KaVo TWAIN Interface                                                     | Romexis                                             |
| KaVo ProXam iCam 1,3m (Mat.no. 2.023.8000)                                                       | KaVo ProXam iCam 1,3m (Mat.no. 2.023.8000) | KaVo ProXam iCam 2,5m (Mat.no. 2.023.8100)                               | KaVo ProXam iCam 2,5m (Mat.no. 2.023.8100)          |
| Operation via camera buttons, treatment unit foot control or touch display of the treatment unit |                                            | Operation via camera buttons, PC or optional Foot Control Vision Full HD |                                                     |
| Optional: connected to PMS software via CONNECTlink (VDDS)                                       | Connected to PMS software (VDDS)           | Connected to imaging software (TWAIN standard)                           | Directly integrated into Romexis (imaging software) |

# KaVo ProXam iCam

## Technical specifications

| Technical specifications KaVo ProXam iCam |                      |
|-------------------------------------------|----------------------|
| Image resolution                          | 800 x 600 pixel      |
| Focus range                               | 9 – 35 mm            |
| Light source                              | 2 white LEDs         |
| Cable length                              | 2.5 m / 1.3 m        |
| Dimensions                                | 176 x 22.5 x 21.7 mm |
| Weight                                    | 33 g                 |
| Power supply                              | 5 V                  |

Thank you for your attention!